


# HI719

## Magnesium Hardness

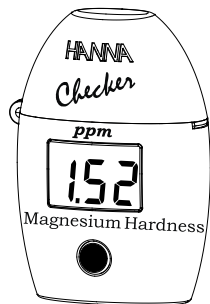

Dear Customer,

Thank you for choosing a Hanna Instruments Product.

Please read this instruction manual carefully before using the instrument. If you need additional technical information, do not hesitate to e-mail us at [tech@hannainst.com](mailto:tech@hannainst.com) or view our worldwide contact list at [www.hannainst.com](http://www.hannainst.com).

### Preliminary examination

Please examine this product carefully. Make sure that the instrument is not damaged. If any damage occurred during shipment, please contact your local Hanna Instruments Office. Each HI719 meter is supplied complete with:

- Sample Cuvettes and Caps (2 pcs.)
- Reagents set for 50 tests
- 1.5V AAA Battery (1 pc.)
- Instruction Manual and Quick Reference Guide




**For more details about spare parts and accessories see "Accessories".**

### Technical specifications

Range	0.00 to 2.00 ppm
Resolution	0.01 ppm
Accuracy	±0.20 ppm ±5% of reading @25 °C/77 °F
Light Source	Light Emitting Diode @525 nm
Light Detector	Silicon Photocell
Method	Adaptation of the <i>Standard Methods for the Examination of Water and Wastewater, 18<sup>th</sup> edition</i> , EDTA colorimetric method. The reaction between magnesium and reagents causes a reddish-violet tint in the sample.
Environment	0 to 50 °C (32 to 122 °F); max 95% RH non-condensing
Battery Type	1.5V AAA (1 pc.)
Auto-Shut off	After 10 minutes of non-use
Dimensions	86.0 x 61.0 x 37.5 mm (3.4 x 2.4 x 1.5")
Weight	52 g (1.84 oz.)

### Functional description

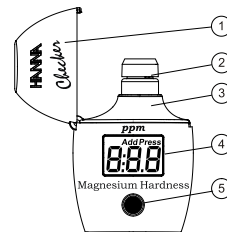

1. Dust cover.
2. Cuvette with cap.
3. Cuvette holder.
4. Liquid Crystal Display.
5. Button.

### Errors and warnings



**Light High:** There is too much light to perform a measurement. Please check the preparation of the zero cuvette.



**Light Low:** There is not enough light to perform a measurement. Please check the preparation of the zero cuvette.



**Inverted Cuvettes:** The sample and the zero cuvette are inverted.



**Under Range:** A blinking "0.00" indicates that the sample absorbs less light than the zero reference. Check the procedure and make sure you use the same cuvette for reference (zero) and measurement.



**Over Range:** A flashing value of the maximum concentration indicates the reading is over range. Dilute the sample and re-run the test.



**Battery Low:** The battery must be replaced soon.




**Dead Battery:** This indicates that the battery is dead and must be replaced. Once this indication is displayed, normal operation of the instrument will be interrupted. Change the battery and restart the meter.



## Measurement procedure

- Turn the meter on by pressing the button. All segments will be displayed. When the display shows "Add", "C.1" with "Press" blinking, the meter is ready.
- Rinse a graduated beaker several times with unreacted sample, before filling it to the 50 mL mark with sample.
- Add 0.5 mL of HI719A-0 reagent, swirl to mix.
- Add 0.5 mL of HI719B-0 reagent, swirl to mix.
- Use this solution to rinse 2 cuvettes. Fill both cuvettes up to the 10 mL mark.
- Add 1 drop of HI719C-0 reagent to one cuvette (#1), replace the cap and invert the cuvette several times to mix. This is the blank.
- Add 1 drop of HI719D-0 reagent to the second cuvette (#2), replace the cap and invert the cuvette several times to mix. This is the sample.
- Place the blank (#1) into the meter and close the meter's cap.
- Press the button. When the display shows "Add", "C.2" with "Press" blinking the meter is zeroed.
- Remove the cuvette.
- Insert the second cuvette (#2) into the meter and press the button.
- The meter directly displays the concentration of magnesium hardness ( $\text{CaCO}_3$ ) in ppm. The meter automatically turns off after 10 minutes.

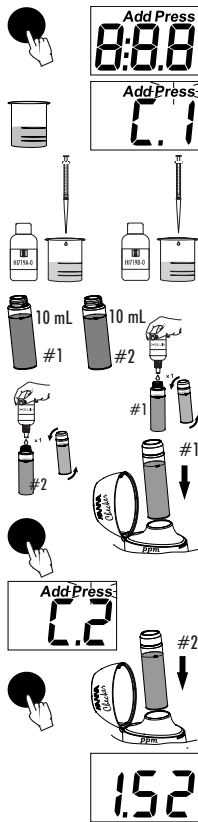

## Tips for an accurate measurement

- It is important that the sample does not contain any debris.
- Whenever the cuvette is placed into the measurement cell, it must be dry outside, and completely free of fingerprints, oil and dirt. Wipe it thoroughly with HI731318 or a lint-free cloth prior to insertion.
- Shaking the cuvette can generate bubbles, causing higher readings. To obtain accurate measurements, remove bubbles by swirling or by gently tapping the cuvette.
- Do not let the reacted sample stand for too long after reagent is added, as accuracy will be affected.
- After the reading it is important to immediately discard the sample, otherwise the glass might become permanently stained.


## Battery management

To save the battery, the instrument shuts down after 10 minutes of non-use. One fresh battery lasts for a minimum of 5000 measurements. When the battery is dead the instrument will display "bAd" then "bAt" for 1 second and then turns off.

To restart the instrument, the battery must be replaced with a new one.

To replace the instrument's battery:

- Turn the instrument off by holding the button until the meter shuts off.
- Turn the instrument upside down and remove the battery cover with a screwdriver.

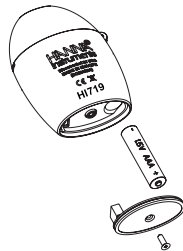

- Remove the battery from its location and replace it with a new one, inserting the negative end first.
- Insert the battery cover and replace the screw with a screwdriver.

## Recommendations for Users

Before using this product, make sure that it is entirely suitable for your specific application and for the environment in which it is used.

Operation of this instrument may cause unacceptable interferences to other electronic equipment, thus requiring the operator to take all necessary steps to correct interferences.

Any variation introduced by the user to the supplied equipment may degrade the instrument's EMC performance.

To avoid damage or burns, do not put the instrument in microwave oven. For the safety of you and the instrument do not use or store the instrument in hazardous environments.

## Accessories

### REAGENT SETS

HI719-25 Reagent set for 50 Magnesium Hardness Tests

### OTHER ACCESSORIES

HI719-11 Magnesium Hardness Certified Standard Kit  
HI731225 Cuvette Black Cap for Checker® HC Colorimeters (4 pcs.)  
HI731318 Cloth for wiping cuvettes (4 pcs.)  
HI731321 Glass Cuvettes and Seal Cap for Checker® HC Colorimeters (4 pcs.)  
HI731353 Cuvette seal for checker HC (4 pcs.)  
HI740028P 1.5V AAA batteries (12 pcs.)  
HI93703-50 Cuvette cleaning solution (230 mL)

## Warranty

HI719 is warranted for a period of one year after date of purchase against defects in workmanship and materials when used for their intended purpose and maintained according to instructions. This warranty is limited to repair or replacement free of charge. Damage due to accidents, misuse, tampering or lack of prescribed maintenance is not covered. If service is required, contact your local Hanna Instruments Office. If under warranty, report the model number, date of purchase, serial number and the nature of the problem. If the repair is not covered by the warranty, you will be notified of the charges incurred. If the instrument is to be returned to Hanna Instruments, first obtain a Return Goods Authorization number from the Technical Service department and then send it with shipping costs prepaid. When shipping any instrument, make sure it is properly packed for complete protection.

*Hanna Instruments reserves the right to modify the design, construction or appearance of its products without advance notice.*